




LIFTS

BRUNO VSL-600 *CURB-SIDER*

*Now you can lift and store your fully assembled scooter
in the rear of your van or minivan the easy way*

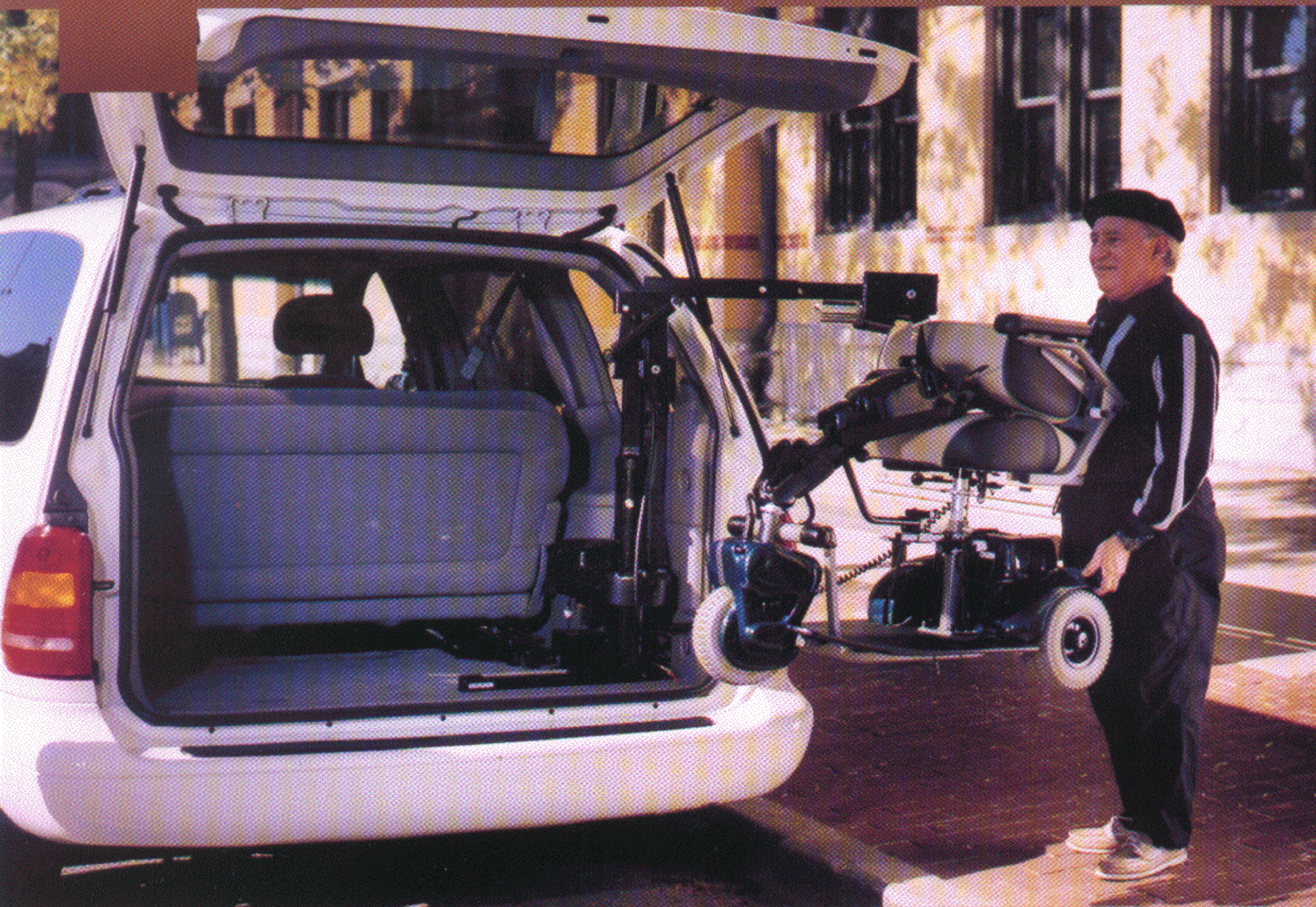

■ The Curb-Sider lift is electrically powered from the vehicle's battery. Simply press a button on the hand-held controller and the Curb-Sider raises and lowers any scooter weighing up to 200 lbs., effortlessly.

■ Bruno's Curb-Sider provides a unique 180° powered rotation, which makes it possible to raise, lower and rotate a fully assembled scooter on the curb at the side of your vehicle, or directly behind it. The Curb-Sider's rotation feature makes it possible to stow scooters in vehicles with tailgates.

■ Bruno's Curb-Sider can end the struggle of manually trying to lift and store a scooter in your van or minivan. The Curb-Sider gives you the convenience you need to transport your scooter and the freedom to travel when and where you like.

■ The Bruno Curb-Sider was designed to fit most types of

vans, minivans, sport utility vehicles and pickups with caps. The Curb-Sider can be adapted for use with some types of power wheelchairs, increasing its versatility.

■ The lift arm can be folded vertically when not in use and installation does not require any structural modifications to your vehicle.

■ The Curb-Sider is only one of 18 different wheelchair and scooter lifts available from Bruno. All Bruno lifts qualify as adaptive equipment and are eligible for reimbursement under the terms of the Mobility Programs offered by GM, Ford and Chrysler.

■ Expert installation and complete after sale service is provided by Bruno's nationwide network of authorized dealers.

Bruno Independent Living Aids, Inc., 1780 Executive Drive, P.O. Box 84, Oconomowoc, WI 53066

All illustrations and specifications in this brochure are based on the latest product information available at the time of publication. Bruno Independent Living Aids, Inc. reserves the right to make changes at any time without notice. FORM LB R VSL-600 1096-2 © Bruno Independent Living Aids, Inc. 1996


Using the Bruno CURB-SIDER is easy!



Just connect the C-arm to the T-bar on the scooter and the Curb-Sider can pick up the scooter from the rear or at the side of the vehicle.



To raise a fully assembled scooter, simply push the lift belt button up.

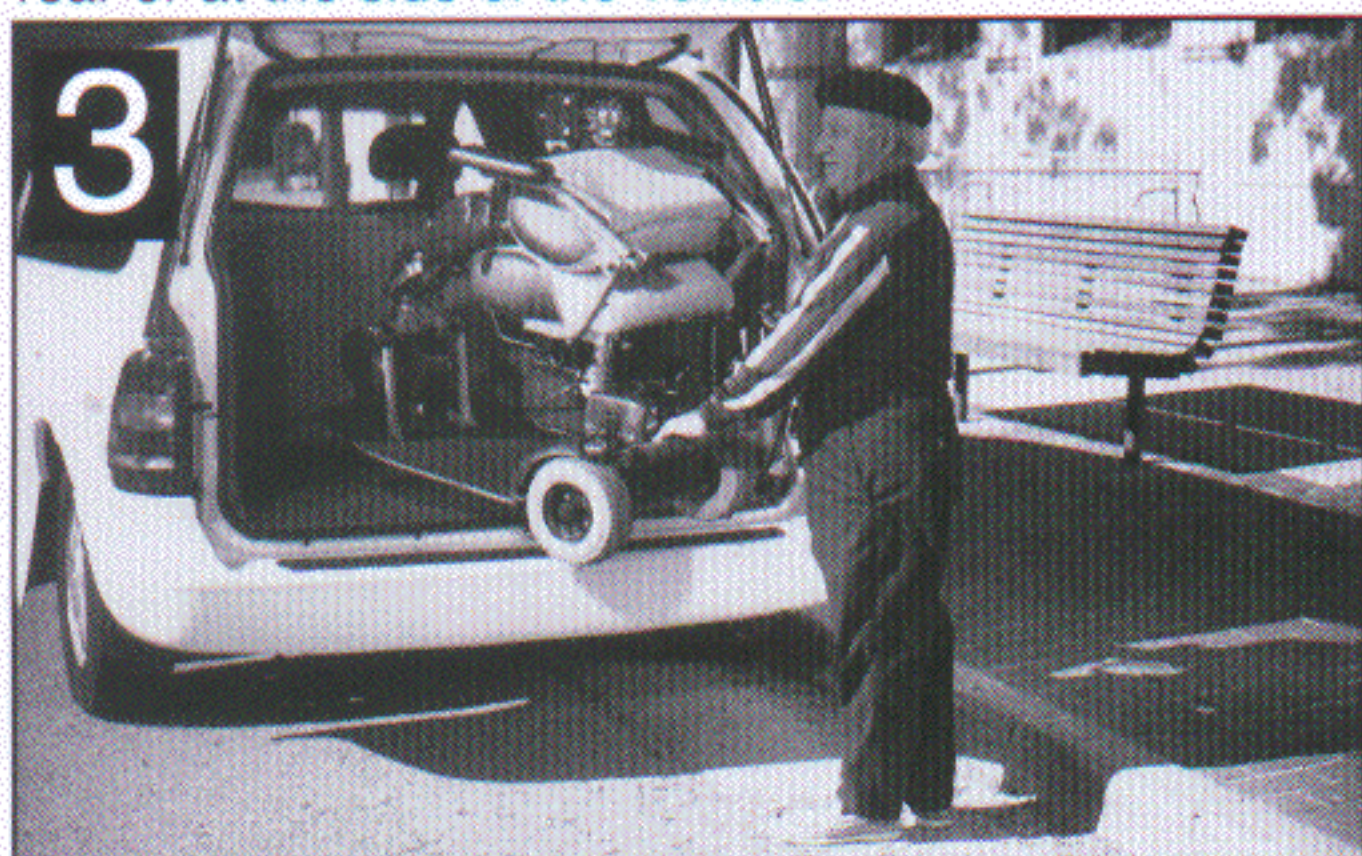


Push the swing in button and the Curb-Sider's motor moves the scooter inward at a comfortable speed.



Simply guide the scooter into the rear opening of your van and position it for storage.



Once the scooter is inside, lower it to the floor. We recommend that you leave the C-arm attached to secure the scooter during travel.

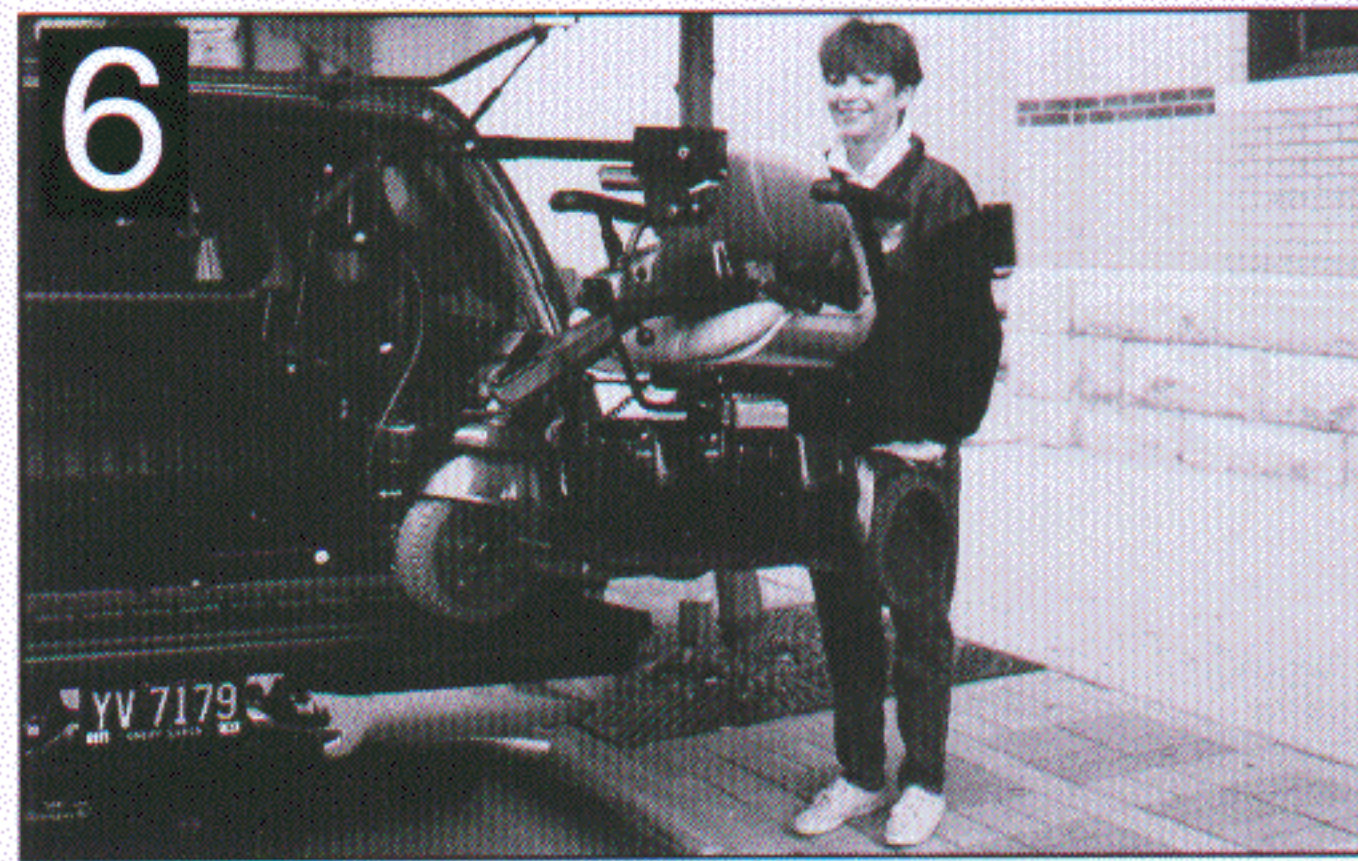

The Curb-Sider can also be used to raise and store a fully assembled scooter in vehicles with tailgates, so long as the opening is large enough to accommodate the lift and scooter.

Represented By:



Call 1-800-882-8183, 1-414-567-4990,
Fax 1-414-567-4341, or <http://www.bruno.com/>